
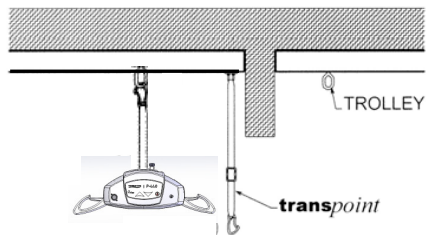


transpoint™

Application Opportunities


The *transpoint* option provides an opportunity to expand the variety of track applications and transfers by allowing a *transportable* lift to pass freely from one track to the next, even under a door header. This gives both an economic and aesthetic advantage where the reconstruction of the door header might otherwise be necessary.

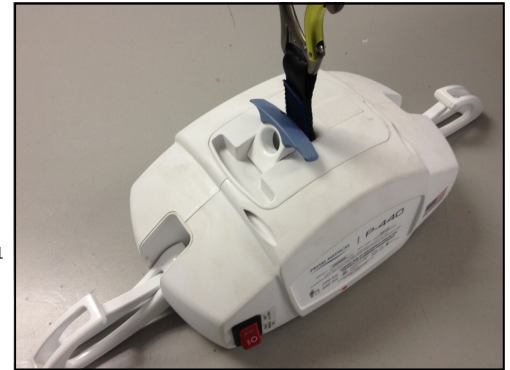

Door Header


Institutions, group homes and other public buildings may require that firewall separations be maintained between all rooms and hallways. With the *transpoint* option, track transfers can be accomplished without modifying door headers.

In bathrooms, the bulkheads separating the shower and/or bathtub areas provide barriers that can be accommodated by the *transpoint* option.

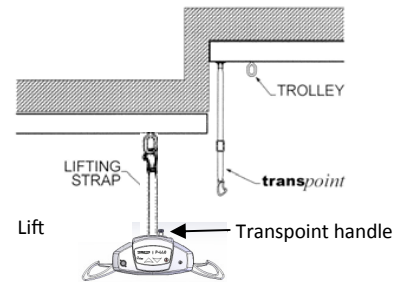

Elevation Change

The *transpoint* can also be used for transfers to tracks at different ceiling heights. The height adjustable *transpoint* strap can be easily modified to fit any combination of ceiling height change, and caregiver height and reach.

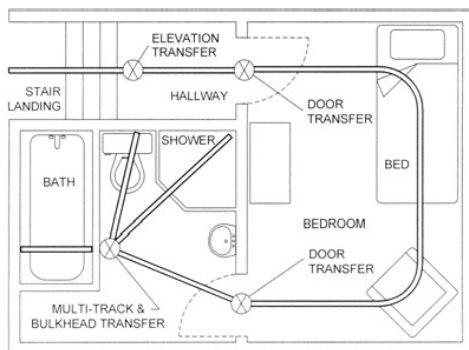

Application Opportunities

Requirements


The *transpoint* option is specified at time of order and consists of two components. Each *portable* (P-200 hoist system is fitted with a handle that can be used to facilitate the transfer hookup. In addition, one *transpoint* transfer strap (adjustable strap and hook assembly) is required for each track transfer point.

The *transpoint* transfer strap comes complete with an end stop that makes the *transpoint* always available for use in multi-track transfers.

Prism Medical UK, Unit 4 Jubilee Business Park, Jubilee Way, Grange Moor, Wakefield, WF4 4TD

Tel: 0844 980 2260 Email: info@prismmedical.co.uk Web: www.prismmedical.co.uk

transpoint – User's Guide to Transfers

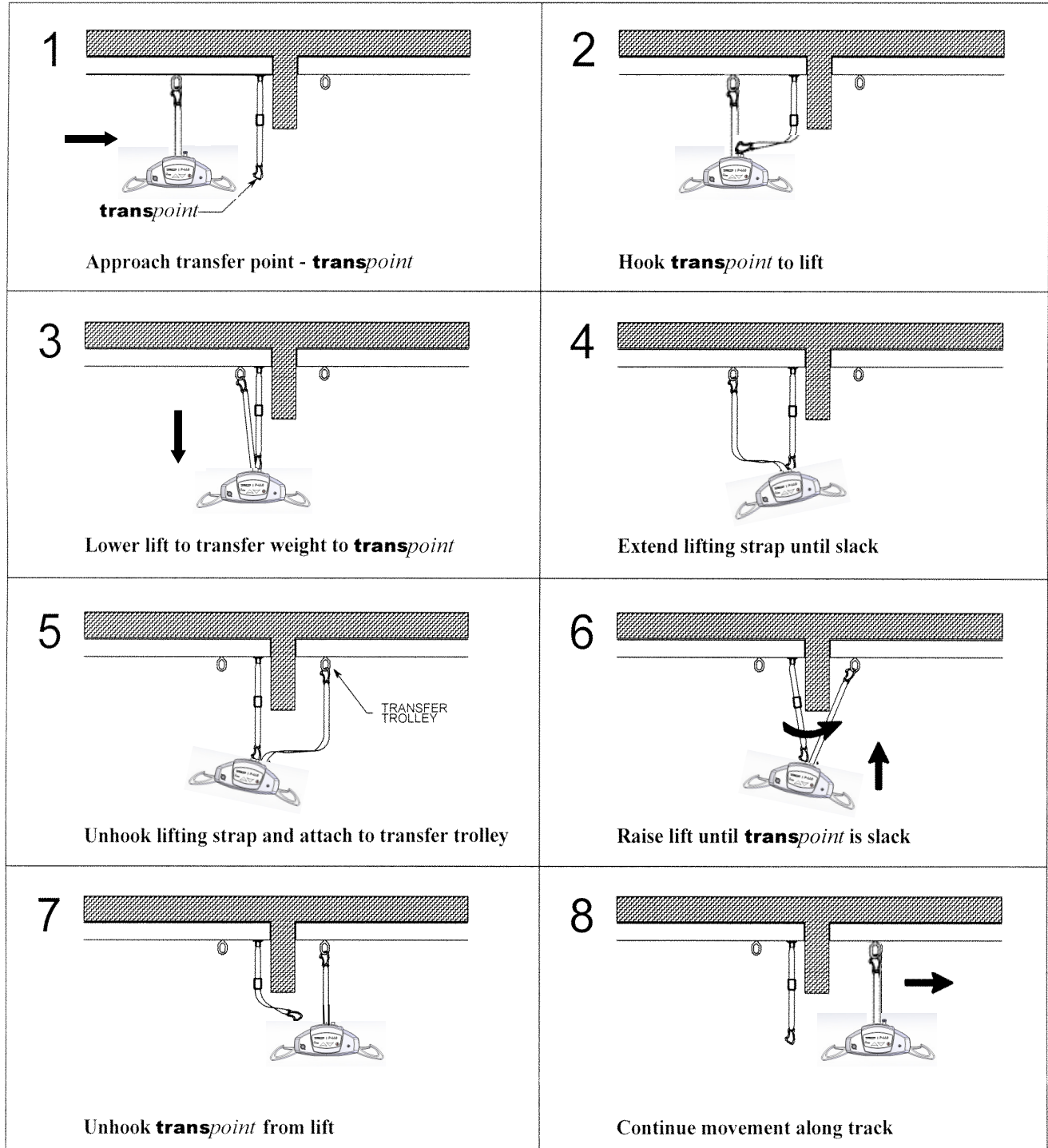

Prism Medical UK, Unit 4 Jubilee Business Park, Jubilee Way, Grange Moor, Wakefield, WF4 4TD

Tel: 0844 980 2260 Email: info@prismmedical.co.uk Web: www.prismmedical.co.uk