
WALK IN BATHTUB



linea
Oceano



OUR RANGE OF "MADE TO MEASURE" PRODUCTS:

The range of Linea Oceano walk in bath tubs for elderly and disabled is really wide:
6 models in 29 different versions to offer the best solution based on the difficulty of movement of the client and the available space in the bathroom.



vasca **Itaca**

the most economical
bathtub



vasca **Tonga**

the ideal bathtub
for small bathrooms

vasca **Samoa**

the bathtub most suitable
for people with disabilities



vasca **Bali**

the bathtub for those
who seek design



vasca **Capri**

the bathtub for those who
want maximum relaxation



vasca **Antigua**


the bathtub for those
who want the top of our range



MATERIAL WE USE:

All our bathtubs are made of hand laid fiberglass and have a shining gel coat surface that is an integral part of the stratification, the supporting frame is made of stainless steel (Aisi304).

CHARACTERISTICS OF OUR BATHTUBS:

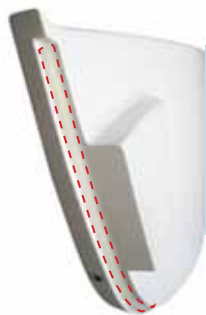

01. SEALING OF DOOR

the door seal perfectly by means of silicone rubber gasket, the main advantage of this material is that the gasket will return to its original shape even after being closed for days greatly reducing wear and tear. The door is fitted with a self-centering locking system.



02. WIDTH OF THE DOOR

can range from 70 to 85 cm, according to the model, obtaining a large, easy and safe entry and exit from the tub.




03. DOOR OPENING

the door opens strictly outwards to 180 ° thus leaving completely exposed the comfortable seating inside (a door that opens inward is a major barrier for those who must get in out of the tub). The door lever is not situated on the inside of the door so there no levers or metallic mechanism of any kind inside the tub that could pose a serious risk of injury to anyone bathing in the tub.

04. INSIDE SEATING

fully exposed built-in seat, wide and 53 cm high. (Not an unstable stool!). It allows you to sit in the tub from the outside; the users rotate their legs inward, like getting in or out of a car. The height makes getting up a minimum effort, as in a chair at home; disabled people can go directly, without gradients, from the wheelchair to the seat inside the tub.

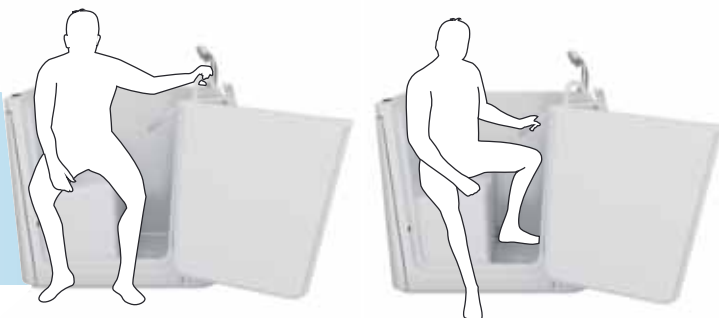

05. ANGLE OF THE SEATING

the particular inclination prevents slipping forward resulting in possible danger (stopping the undesirable slippage caused by the inevitable floating of the leg due to the gradual filling of the tub).



06. RECESS FOR HYGIENE

positioned in the middle of the seat allows the personal hygiene without having to get up from a comfortable position.



07. GRAB BAR

The grab bars placed on board and inside of the tub, gives the user safe and free movement getting in or out of the tub.



08. TAPWARE

located on the front of the tub comfortably near your fingertips, consisting of mechanical mixer, tub overflow, water in let, pull-out hand shower, on demand shower rod.





Via B. Luini, 5 – 20825 Barlassina (MB) Italy
Phone: +39 0362 556354 – Fax: +39 0362 557577
Cell: +39 328 6070195 – Email: marcello@lineaoceano.it

www.lineaoceano.it

